



TENSMED S84



The top
professional product
for pain relief and muscle stim



PARTNER FOR LIFE

TENSMED S84

The top professional product for pain relief and muscle tim

The TENSMED S84 is the ultimate choice for professional treatment of pain related disorders and muscle stimulation in physiotherapy, rehabilitation, sport, fitness and beauty & aesthetics. With this powerful device you can stimulate up to four large muscle groups or different pain areas at the same time. Or give bilateral treatments simultaneously (for instance left and right m. Quadriceps). The most advanced technology is used to guarantee the most user-friendly and safest use of the stimulator in any treatment situation.

With the TENSMED S84 you can treat chronic pain, patients recovering from a traumatic condition, immobilisation as well as maintain or improve muscle function in healthy individuals using TENS (Transcutaneous Electrical Nerve Stimulation) and NMES (Neuro Muscular Electrostimulation) currents. This makes the device ideal for sport and rehabilitation.

The TENSMED S84 is a 4-channel device with a clearly legible LCD display and easy to operate.

265 Programs

- 13 Rehabilitation
- 3 Incontinence
- 33 Pain (TENS)
- 53 Sport
- 29 Fitness
- 16 Beauty & Aesthetics
- 28 Action Now
- 27 FES Now ('3S')
- 48 NMES Now
- 15 Free programs



TENSMED S84

PROGRAMS

Programs

250 Specific preset programs for functional rehabilitation, pain management, sports, fitness and beauty & aesthetics.



Easy program menu

Rapidly create new TENS and NMES programs or modify existing ones using a great variety of stimulation parameters. The stored 'Easy programs' are accessible for all users.



Syncro-stim

Increase/decrease the intensity of all 4 channels simultaneously using only one button.




Run time function

Option to adjust treatment time, frequency and pulse width during the treatment.



My Trainer

Displays images for the correct placement of the electrodes.

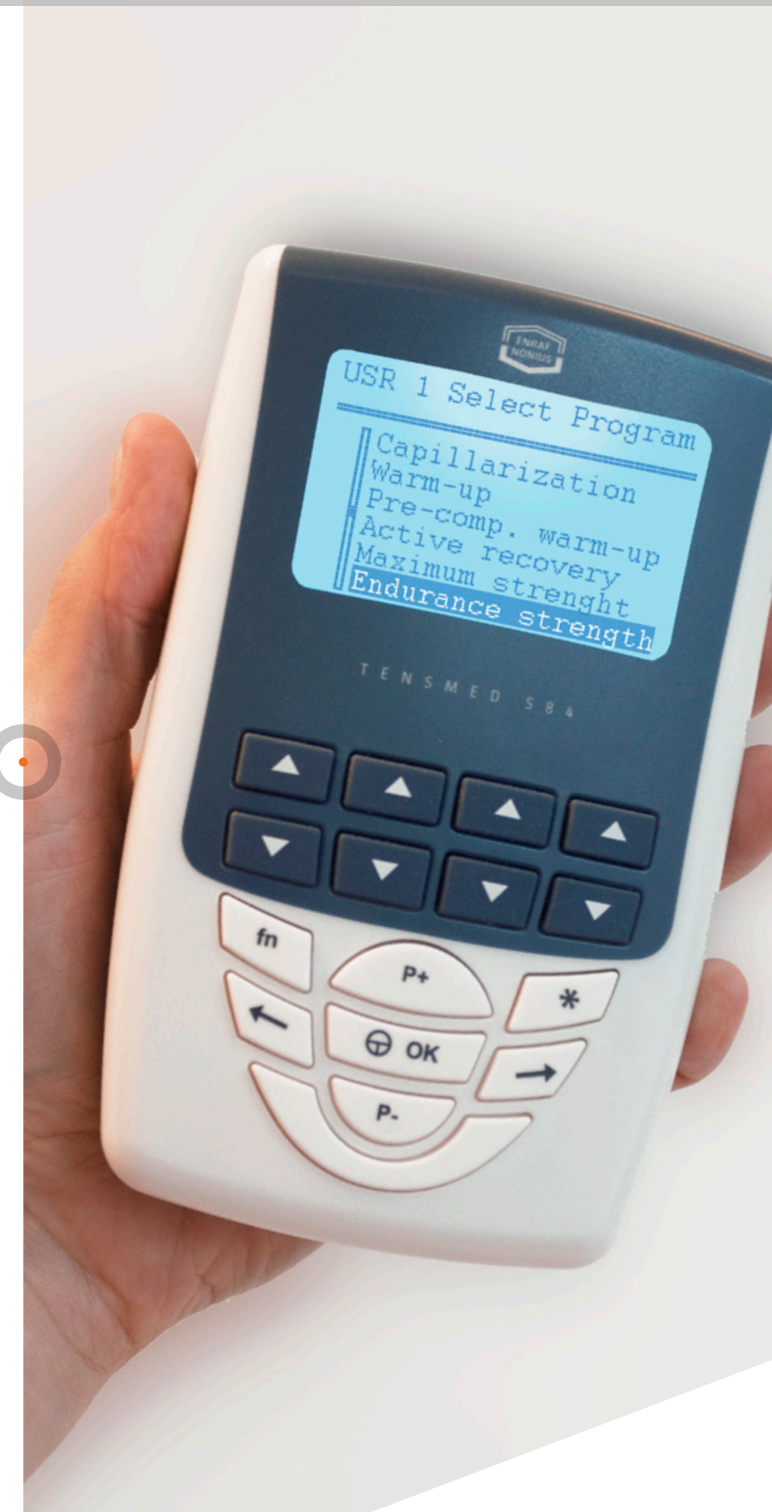



TENSMED S84



Programmable

15 Free memories to create and store your personalized TENS and NMES programs.



Quick menu

Fast access and easy navigation to functional Neuromuscular Stimulation programs: 'Action Now', 'FES Now' and 'NMES Now'.



Action now

Special function that allows you to start a single muscular contraction by pressing a button on the device. Action Now programs are particularly useful to combine and synchronize the electric stimulation with a voluntary action to obtain a greater recruitment of the muscular fibers and to improve coordination.



FES now ("3S")

Special function for 'Serial Sequential Stimulation'. FES Now programs are Functional Electro-Stimulation programs characterized by a delayed activation of channels 3 and 4 compared to channels 1 and 2. They permit the stimulation of the musculature in a kinetic chain. In the aesthetic field, 'FES Now' programs are used for drainage.

TENSMED S84



NMES now

Fast access to specific muscle stim programs with preset ON/OFF stimulation times (intermittent stimulation) .



Multi-user

Up to 10 different users can create and store a library of personalized programs in the 'Favorites' and 'Last-10' memory.




Last 10

Automatically stores the last 10 executed programs. 'Last 10' programs will remain available for a rapid and easy execution.



Favorites

Permits to store up to 15 customized TENS and NMES programs for each user.

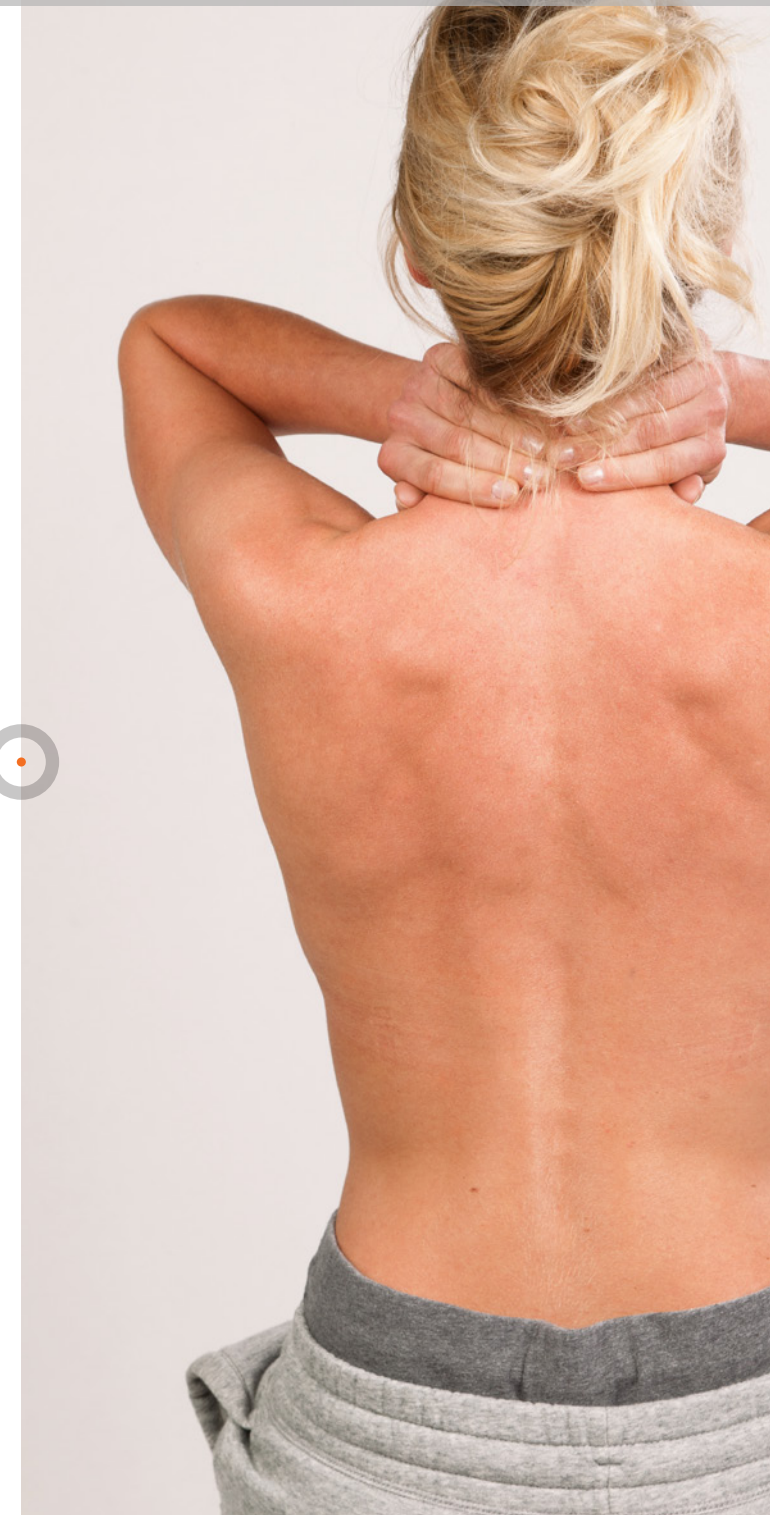



TENSMED S84



Auto-stim

Execute a program in automatic mode (= without having to regulate the intensity manually). Intensity values are automatically set to the levels used during the last execution of programs listed among the "Last 10" memory.



2 + 2 mode

Allows you to execute two different programs (NMES or TENS) at the same time, thus permitting the simultaneous treatment of two patients or two muscular groups.



Working time

Indicates the total time the device has been used for stimulation treatment.



Backlit display

The clearly legible LCD display contributes to optimal usability and ease-of-operation.



TENSMED S84 programs

Rehabilitation

Muscle stimulation programs are particularly indicated for the functional recovery after surgical articulation interventions and after immobilization time. Programs included for Vastus Medialis Reinforcement and Recovery after ACL Surgery, Prevention Programs (prevent traumas and injuries), Re-Athletization and Hemiplegia.

Incontinence

Specific programs for the treatment of urological disorders such as stress, urge and mixed incontinence are included in the Rehabilitation section.

Pain

Eliminate pain with the special pulses of the TENS currents that, stimulating your body in a soft way, enables you to recover your health, wellness and a renovated life enjoyment. With the TENS programs of the TENSMED S84 you can treat acute and chronic pain that is related to musculoskeletal disorders. Programs included for Acute and Chronic Pain, Cervical pain, Lumbago, Periarthritis, Epicondylitis, Sciatica, soft tissue injuries, etc.

Sport

The TENSMED S84 is programmed in order to selectively stimulate fast and slow fibers to train different muscle qualities. Includes programs for Warm-up, Maximum strength, Endurance strength, Explosive strength, De-contracting (muscle relaxation), Capillarization (increases the blood flow to the muscle and improves endurance and recovery abilities) and Pre- & Post-competition/ training and Active Recovery (reduce muscle fatigue after training or competition).

TENSMED S84 programs

Fitness, Beauty & Easthetics

With the TENSMED S84 you can care for your body and defend its vitality and beauty. Use the specific beauty programs to treat skin flaws: improvement of the capillary microcirculation, toning and tissue firming, liquid drainage, increase in metabolism and cellular activities, reduction and mobilization of the fat in the localized deposits (e.g. cellulite, water retention, localized fat).

The TENSMED S84 is also ideal for Loss of weight/ Localized lipolysis, Anti-cellulite treatment, Lymphatic drainage and Firming/Toning. The Firming and Toning programs contribute to reduce skin flaws through the general and localized increase in the metabolism.



TENSMED S84 programs

QUICK NMES (Neuromuscular Electrostimulation)


The Quick NMES menu offers 3 ways to stimulate and train muscles (or muscle groups) in the most functional way, using one of the 103 dedicated programs. In just 3 steps you can select one of the following Quick NMES options: 'Action Now', 'FES Now' or 'NMES Now'. Fast, easy and functional.

Action Now programs are normal NMES programs with the only difference that you start the contraction manually: after pressing a button of the TENS MED S84. These programs are especially useful to link and synchronize the electrical stimulation to a voluntary contraction. The programs vary from very rapid and brief contractions (typical for explosive movement) to slower and longer contractions (when aiming at hypertrophy).

FES Now ("3S") programs are characterized by an activation delay of the channels 3 and 4 compared with the channels 1 and 2. These programs permit to stimulate the musculature in a kinetic chain thanks to the differentiated activation times of the muscular groups involved. In the aesthetic field, the 3S programs allow to create a real sequential drainage. The rhythmic stimulation of the musculature causes an increase in the speed of the venous circulation (hematic flow) and lymphatic circulation (lymphatic flow)

NMES Now offers fast access to specific muscle stim programs with a preset ON/OFF stimulation time.

There is an alternation between working time and rest time (active and passive). For example, 6 seconds of contraction and 10 of recovery, During the time of recovery, there is also the possibility to raise the intensity of current in order to execute a cool-down active recovery during the rest.

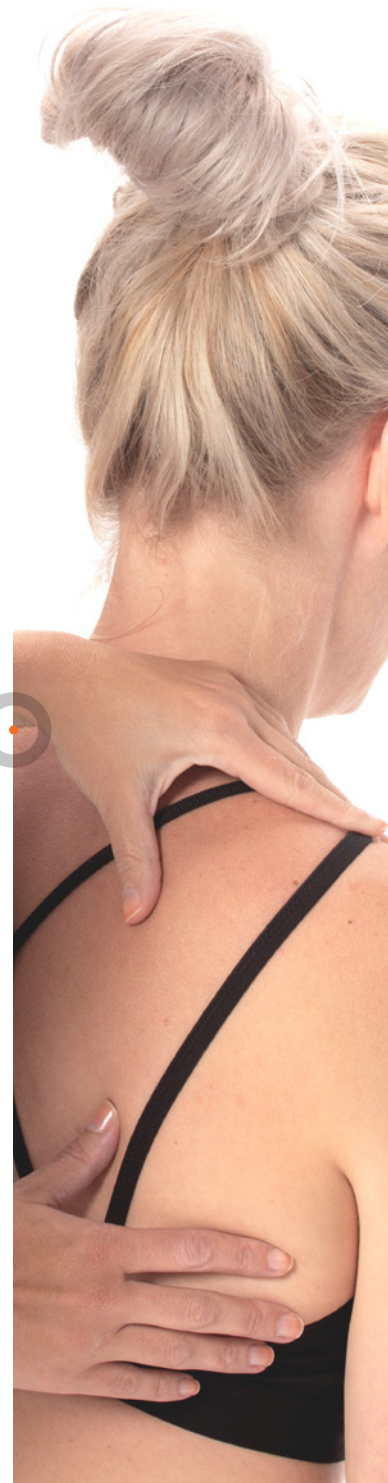

TENSMED S84

1727911 Tensmed S84

Standard accessories

Carrying bag

4 Cables for electrode connection


4 Self-adhesive square electrodes (50x50mm)

4 Self-adhesive rectangular electrodes (50x90mm)

Charger

User manual

ORDERING DATA

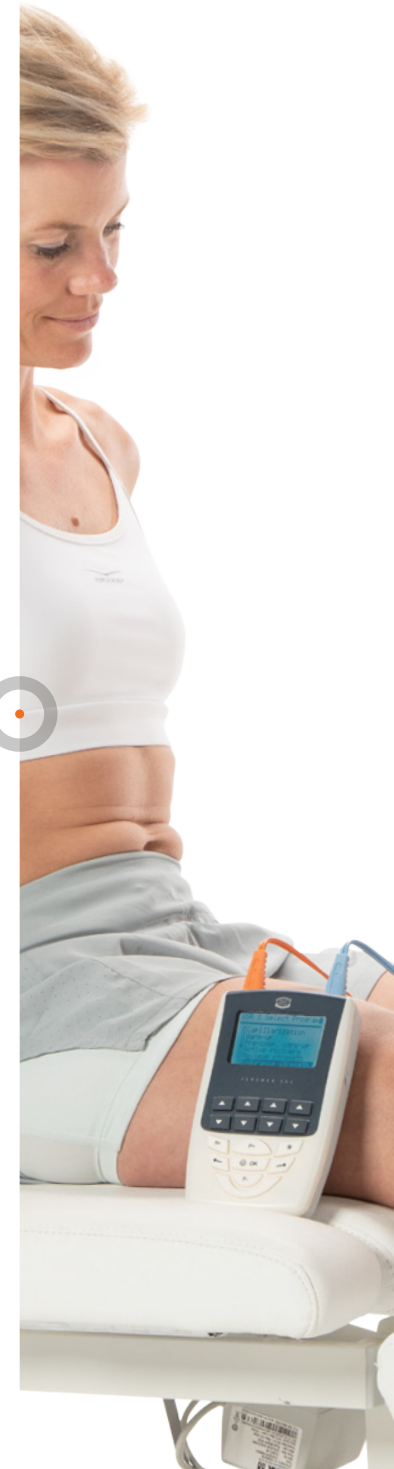


TENSMED S84

TECHNICAL SPECIFICATIONS

Numbers of channels:	4 independent channels
Intensity:	0-120 mA per channel
Wave form:	symmetrical biphasic pulse, 100% compensated
Constant current:	up to a resistance of 1000 ohm
Frequency:	0.3 – 150 Hz
Pulse width (single phase):	50-450 μ s
Total pulse width	100-900 μ s
Number of pre-set programs:	250
Numbers of personal programs:	15
Stimulation forms:	Continuous / conventional, Burst, Frequency Modulation, Pulse Width Modulation Phase Duration Modulation), Constant, Intermittent, Active Rest (electrostimulation during resting phase)
Timer:	max. 120 min (4 subsequent phases of each 30 min)
Power supply:	Battery
Battery:	Rechargeable battery (NI-MH / 1.8 Ah / 7.2V)
Charger:	Input: 100-240VAC~50-60Hz  0.2A Output: +8.8V \square 0.2A
Dimensions:	160x99x35.4 mm
Weight:	404 g
Display:	Backlit display

This is a medical device.

It complies with all applicable medical regulations. For more details please consult the instructions for use (downloadable via www.enraf-nonus.com).

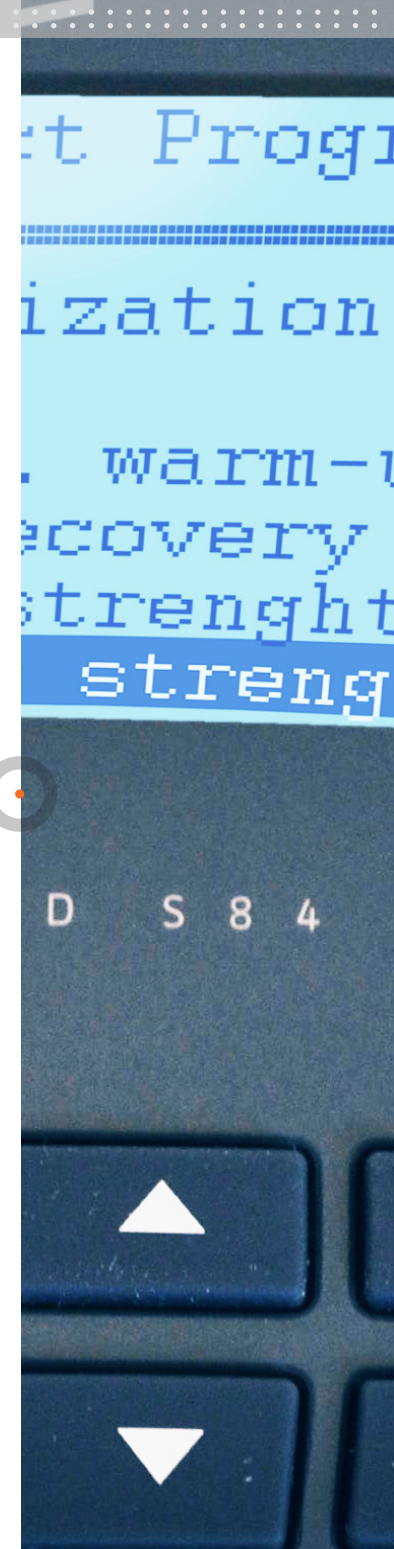


TENSMED S84

ACCESSORIES

Adhesive electrodes


- 3444222 EN-Trode Ø 2,2 cm, 2 mm female, 10 sheets of 8 (also for EMG)
- 3444056 EN-Trode Ø 3,2 cm, 2 mm female, 10 sheets of 4
- 3444135 EN-Trode Ø 5,0 cm, 2 mm female, 10 sheets of 4
- 3444057 EN-Trode 5x5 cm, 2 mm female, 10 sheets of 4
- 3444058 EN-Trode 5x9 cm, 2 mm female, 10 sheets of 4

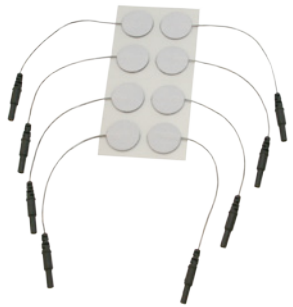


3444222



3444056

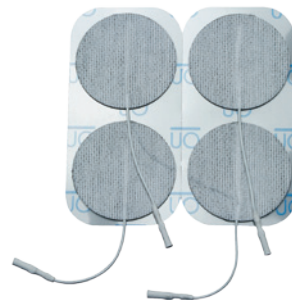


3444135

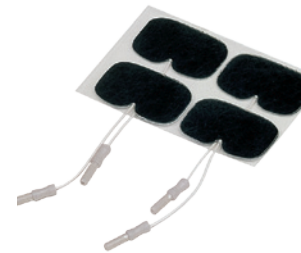

3444057



3444058

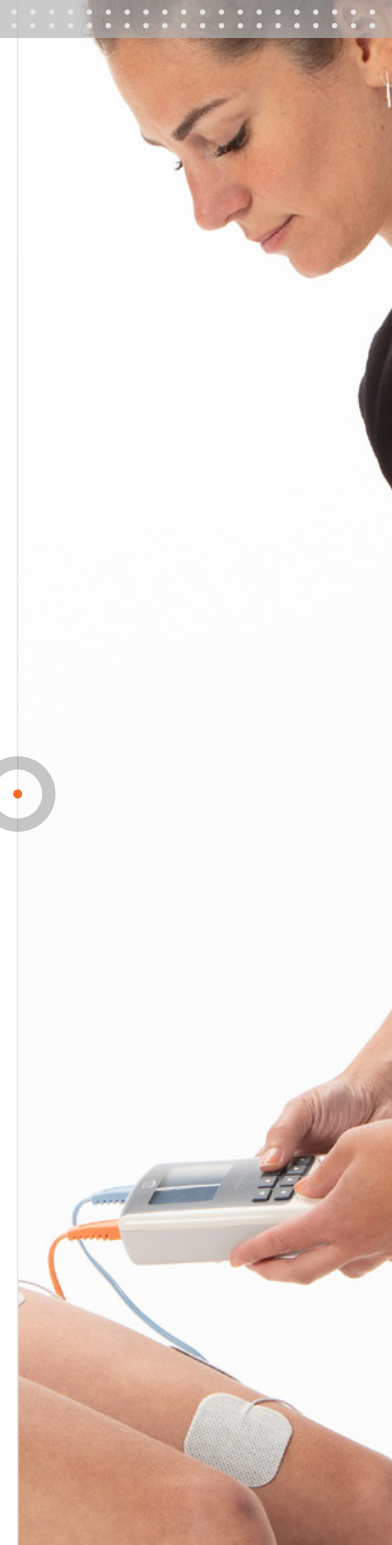

TENSMED S84

ACCESSORIES

Rubber electrodes

- 3444128 Rubber electrodes 4x6 cm, 2 mm female, set of 2.
- 3444129 Rubber electrodes 6x8 cm, 2 mm female, set of 2
- 3444130 Rubber electrodes 8x12 cm, 2 mm female, set of 2
- 3444380 Silicon electrode 5x5 cm TENS, 2 mm female



3444128




3444129



3444130



3444380

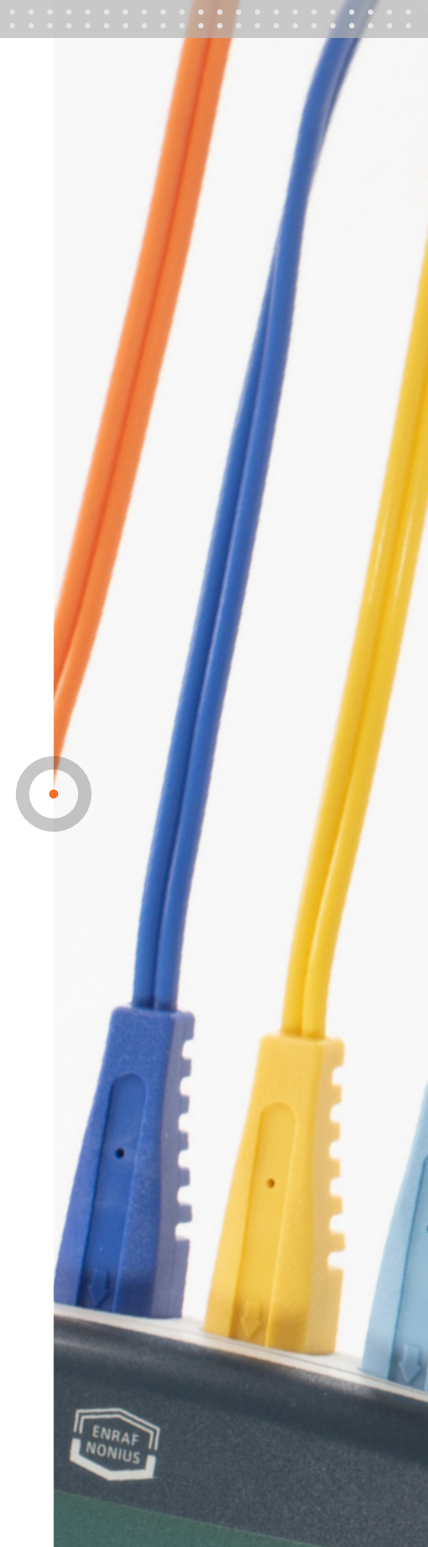

TENSMED S84


ACCESSORIES

Moist pads for rubber electrodes


- 1460273 Moist pads for rubber electrodes 4x6 cm, set of 4
- 1460266 Moist pads for rubber electrodes 6x8 cm, set of 4
- 1460275 Moist pads for rubber electrodes 8x12 cm, set of 4

Adapters


- 2523524 Adapter plug, 2 mm female, 4 mm male, red
- 2523523 Adapter plug, 2 mm female, 4 mm male, black

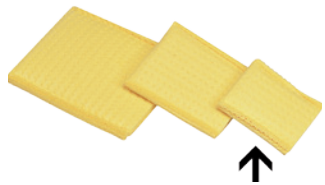

1460273

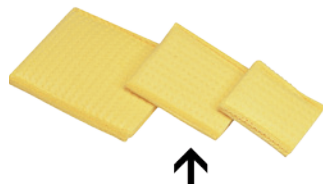

1460266

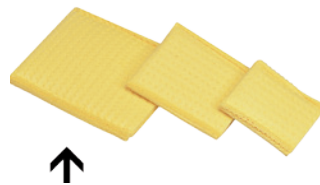


1460275



2523524



2523523

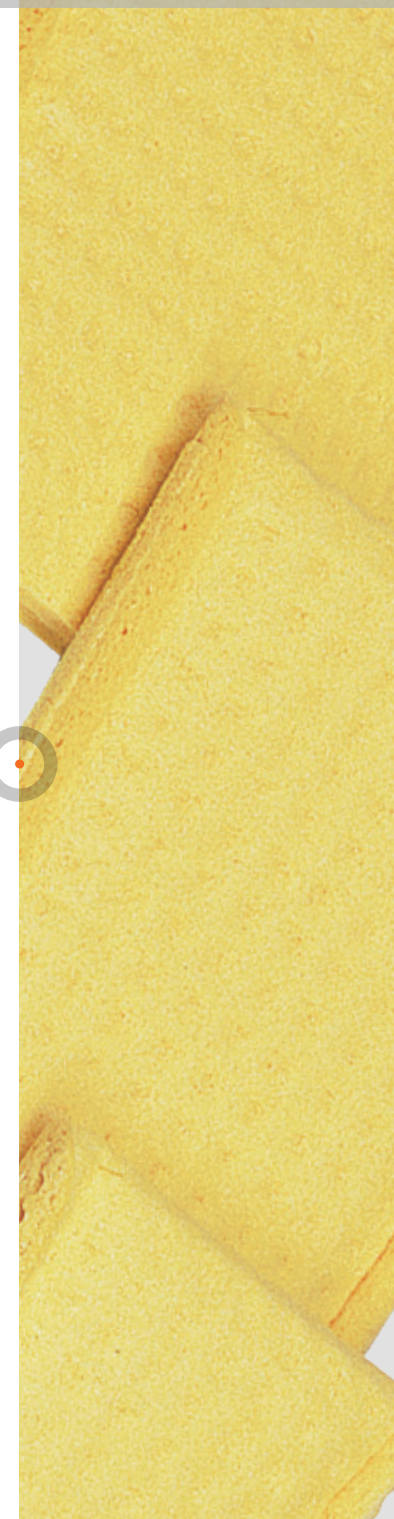



TENSMED S84

The top
professional product
for pain relief and muscle stim